
## KULLANIM KILAVUZU

### Elektronik Akıllı Tartı

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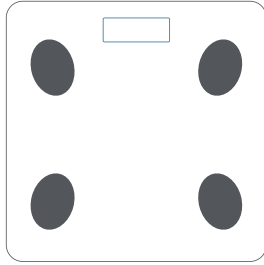[www.simart.me](http://www.simart.me)

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
### ŞİMART Application

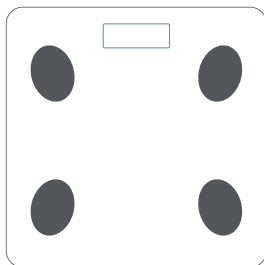[www.simart.me](http://www.simart.me)

**IMPORTANT:**  
Please read the manual carefully before using the device. Failure to comply with the instructions or using this scale in ways other than the ones mentioned in this manual may result in injury to yourself or damage to the scale.

#### PARAMETERS

• Weight	• BMI	• Lean Body Weight
• Body Fat	• Body Water	• Muscle
• Calories	• Visceral Fat	• Protein
• Bone Mass	• Body Type	• Body Age

#### INDICATION SYMBOLS

Lo	Low Battery - Please Change the Batteries
Err	Overload
C	An error has occurred during measurement. Step off the scales and step on them again to repeat the measuring process

#### ! WARNING

- DO NOT use with medical implants such as pacemakers.
- DO NOT stand on the edge of the scale or jump on it.
- DO NOT overload the scale (MAX. 180kg).
- DO NOT drop scale or drop an objects on it as this may damage the sensors.
- DO NOT immerse the scale in water or use chemical cleaning agents. Clean the scale with slightly damp cloth.
- NOT recommended for pregnant women.

#### INSTRUCTION

##### 1. User Registration

- 1) Download the Şımart APP in application stores, or scan the QR code below.

- 2) Create a new account with email address
- 3) Get the verification code to finish the registration.

##### 2. Bluetooth Connection

- 1) Please open the Bluetooth and GPS on phone, click "+Add device" to automatically scan device (Note: please keep the scale running during connection)
- 2) For Wi-Fi version, please click "Configuring Wi-Fi" to connect.
- 3) After connection, click the device to enter the measuring page.

##### 3. Start Measurement

- 1) Select the user and fill in information with the photo profile on up-left corner
- 2) Step to active the scale and check if the data is in sync with the phone.
- 3) Step again with bare feet and keep the body upright and wait for the measurement finished.
- 4) For the first time measurement, please click data to claim on account

#### 4. More Functions

##### 1) Weight Unit Switch

Change the unit between kg/lb/st.

##### 2) Data Trend

Check Weight, Body fat and BMI in chart by week, monthly or yearly data.

##### 5. Battery Advice

- If the scale is not to be used for long periods, it's recommended to remove the battery to avoid damage to the scale due to possible battery leakage.
- Keep the battery out of reach of children.
- Do not mix old and new batteries, with different compositions or of different brands in order to prevent possible leakage, explosion.
- Do not heat or deform the batteries or explore to fire.
- Waste batteries should not be disposed of with household waste.
- Please check with your local authority for battery recycling advice.
- Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice.

#### Disposal

The packaging material is recyclable. Please dispose of it in an environmentally friendly manner and take it to a recycling centre.

Dispose of the device in an environmentally friendly manner. Dispose of it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

The batteries should be removed prior to disposal of the device and disposed of separately from the device. To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but must be taken to suitable collection points. Also note the relevant statutory regulations governing the disposal of batteries.